

Safe and careful patient movement



Venture Company
Yu Gwang Precision



미국의료기기



호주의료기기



의료기기



특허



벤처기업



ISO






단체표준



우수제품



자율안전

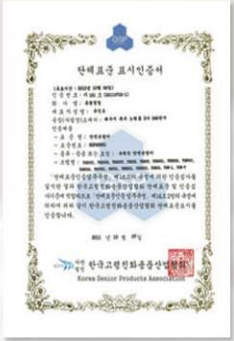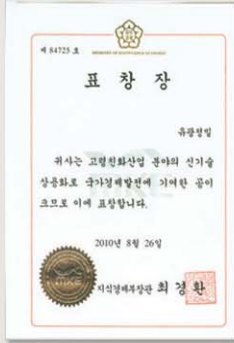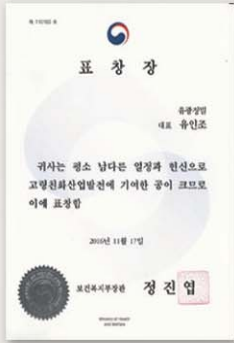



Lift Chair

YGM005



User of Lift Chair

1. Standing Up

For people who are able to walk once they are helped to stand up

2. Sitting

For helping people to sit down on the floor or moving from one seated location to another.

3. Getting on and off a wheelchair

For people who cannot get on or get off a wheelchair by themselves. Height adjustment allows the assistants to easily help the patient.

4. Everyday Lift

Getting on and off a toilet seat, getting close to the height of most indoor furniture, and moving around indoors.

5. Work Place

At the work place, adjusting height while working allows for better posture and comfort.

How to use and operate

1. Secure the lift by pressing the safety device on the rear wheels
2. Draw the armrest completely to the back
3. User moves onto the seat
4. Activate the lift with remote controller S/W pressing up and down.

Notice

- Watch out for your legs when going down—make sure it does not get stuck under the chair (*Switch would not work because of the safety device attached)




Product specifications

Power	AC 110/240V 50/60Hz
Size	800(W) × 600(L) × 1050(H)mm
Weight	50kg
Capacity	120kg
Material	Steel, Polyester

Transfer Lift

YGM004



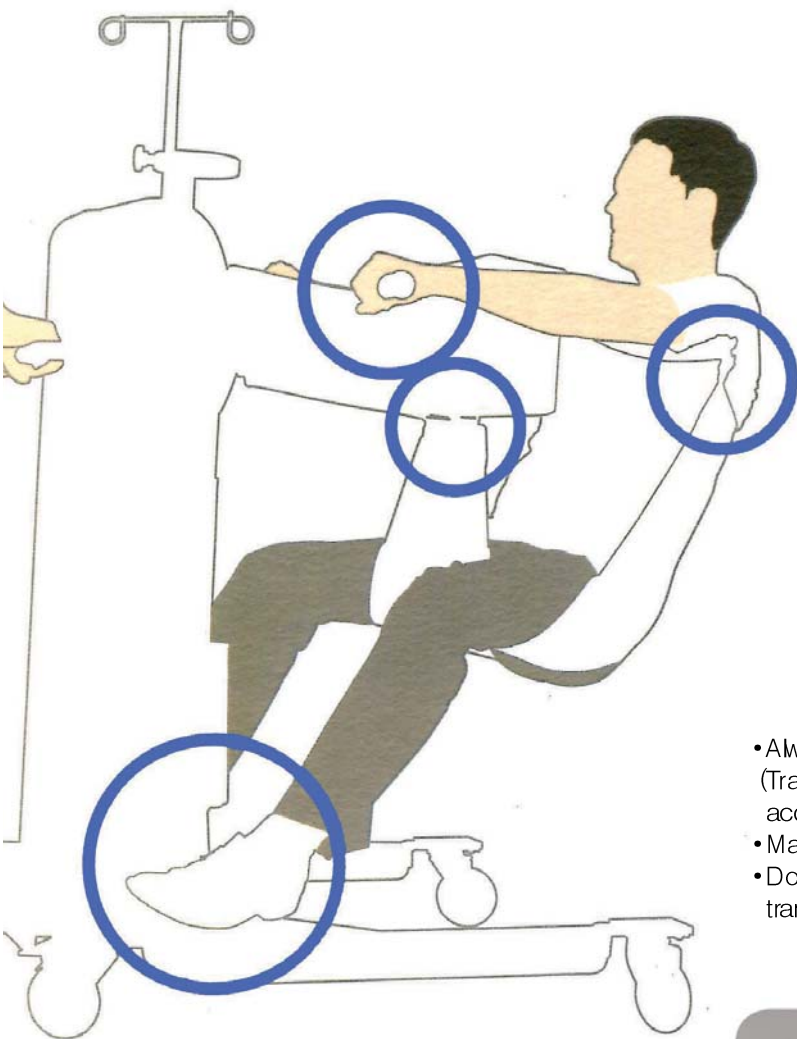

Easy to use
Simple, Convenient steps
One person operation
Lifts up to 400 lbs



- Always hold onto the handle during transfer. (Transferring the patient without holding the handle may cause accidents)
- Make sure the sling hook holders are fastened securely.
- Do not swing legs and place both feet on the footrest during transfer.



Product specifications

Power	AC 110/240V 50/60Hz
Size	1,015(W) × 750(L) × 1,195(H)mm
Weight	74kg
Capacity	180kg
Material	Steel, ABS copolymer, Polyester



design registration




senior-friendly products



Group standard specification


Transitions modular entry ramp

It helps to move and solves getting caught in the threshold people with disabilities who use transfer patient lift, wheelchair, mobility scooter, support walker

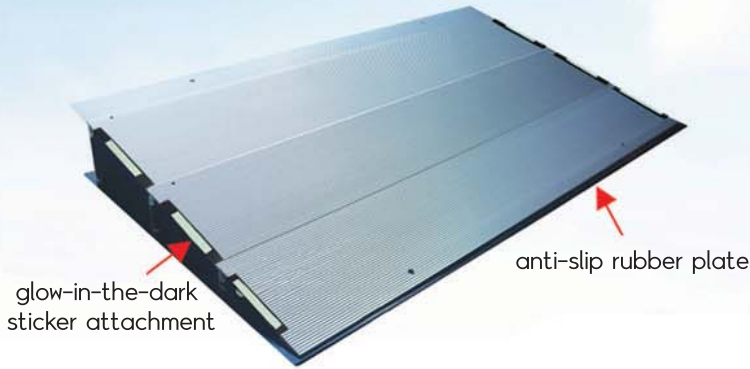


YGM2

Group standard number : KSAP 1019
 verification code : 729(20200309-1)
 senior-friendly products : F24000045001

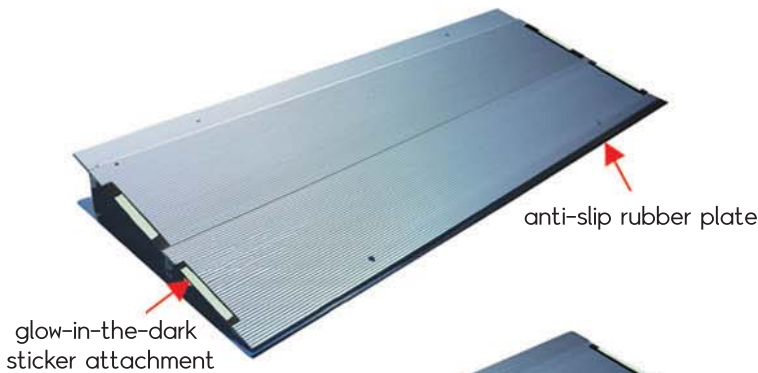


YGM3

Group standard number : KSAP 1019
 verification code : 730(20200309-2)
 senior-friendly products : F24000045002

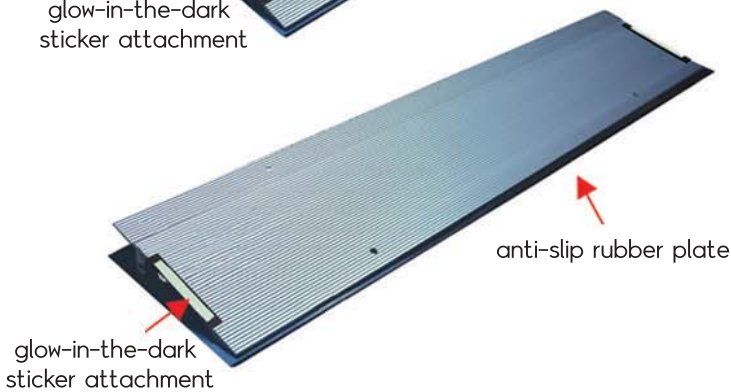

YGM4

Group standard number : KSAP 1019
 verification code : 731(20200309-3)
 senior-friendly products : F24000045003

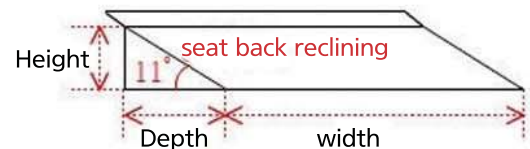

Product specifications

Item	Model	Size	a door sill height	weight	allowable load
		width x depth x height(mm)			
Transitions modular entry ramp	YGM2	750 x 420 x 75	65-81	5.5kg	100kg
	YGM3	750 x 290 x 50	35-56	3.48kg	100kg
	YGM4	750 x 165 x 25	12-31	1.83kg	100kg

Quality of the material : aluminium, acrylonitrile butadiene styrene copolymer

Composition : body, driver, cap, anti-slip rubber plate, bolt, plastic anchor

caution in use

- Before use, please check if there's anything wrong.
- When fixing to the floor, use double-sided tape, bolt, plastic anchor, etc. to secure.
(If the ramp is not fixed to the floor, it may move and cause injury, such as slipping, etc.)
- When constructing the floor, pay attention to damage of water supply, sewage, heating, and pipes.
- Do not use more than 100 kg of safety operation load.
- Assemble and use it according to the height of the place you want to use it.

assembly method

When assembling, gently push into the groove as shown below.



How to Install


- Remove any foreign substances on the floor you want to install.
- Remove the double-sided tape film that is attached to the underside of the slope and attach it to the floor.
- Fix it using wood or concrete bolts, and plastic anchor.
- If construction is difficult due to water leaks, such as toilets, entrances, etc., or damage to pipes, it can be used on a rubber plate.



Securing with M4 bolts to prevent dislocation



Remove and install double-sided sticker film

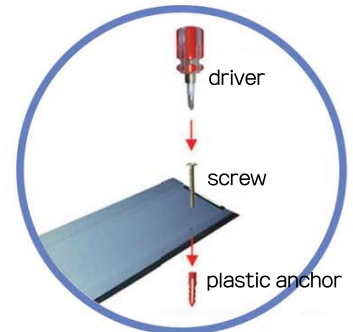

Fix it using bolts, and plastic anchor.

Portable Ramps



Precautions

- Scooter, manual wheelchair and power wheelchair users will gain new access where elevated areas impede mobility.
- Lightweight, durable threshold ramp installs easily and makes entering or exiting convenient and safe.
- Works with sliding doors, thresholds, raised landings, etc.
- Measure door entry height and choose correct ramp—use, either just one or a combination of two depending on your needs.
- Anodized aluminum is rustproof and corrosion proof with a non-skid surface.



Product specifications.

Size	785(W) × 430(L) × 82(H)mm
Capacity	300kg
Material	Aluminium, ABS copolymer
Weight	6.2kg

How to assemble



Examples of use



Darkness



Test



Examples of use

- 1 Before use, please check whether or not each part is working.
- 2 Secure the body using M4 bolts and 3mm pins to prevent breakaway.
- 3 When fixing to the floor, please use double-sided tape, silicone, anchor bolts, etc.
- 4 Safety Operation Load: Do not allow weight to exceed 300kg.
- 5 Assemble according to the height of the place you are using it in.
- 6 To prevent scratches and slips on the floor, please use the rubber plate to assemble and install the ramp.
- 7 When assembling, slide it gently into the groove.

Urinecare¹

(for men)

For people who require assistance when urinating at night.
Used with the "pessmo" cover on penis.



Components



Why Urinecare is beneficial



It's Convenient.

Caregivers are no longer disrupted at night



It's economical.

The need to buy expensive diapers is gone.



It's hygienic.

The urine is stored in the vinyl hose, which allows the urine to be discarded and thus reducing contamination.



It's environmentally friendly.

It reduces the amount of waste from diapers.

How to use Urinecare



*Don't worry.
It does not slip away.*



1 Wrap the "Pessmo" cover around the shaft of the penis.

2 Cover the "Pessmo" cover with PE film, and tie it up with the Gismo band.

3 Place the anti-reflux tube under the knee and pass the PE film through.

4 Fit the PE film into the Urinecare lid hole.



Design Registered



S-line



Pelvis balance



Hips up



Strengthening nether limbs

Keegel Flying Wing

Make yourself healthier and beautiful by investing 10 minutes daily!

Built-in count function measure the time of exercise correctly so that you can not only know your daily amount of exercise but set your target count for exercise conveniently.

Built-in Counter!

- Convenient checking of exercise count
- Setting target by checking daily amount of exercise!

01 Lighter than before!




Upgraded resilience core

Upgraded durability by carbon steel material!!!
Available to use longer
Improved exercise effect!!!


02



Experience "Pelvis correction" now by Keegel Flying Wing~!!

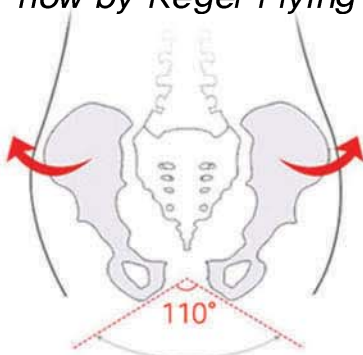

Before Exercise Widened Pelvis

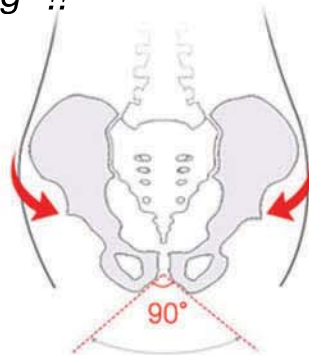

After Exercise Contracted Pelvis

Exercise strengthening Pelvis Muscle

WOMAN

- Preventing damage on perineal region in childbirth
- Fast recovery after childbirth
- Strengthening muscle of perineal region
- Preventing incontinence of urine
- Supporting womb

MAN

- Preventing prostatitis
- Preventing premature ejaculation
- Improving stamina

Product Specification

Product Size	W 20cm X L 7.2cm
Material	SK5, Moplen
Color	Grey

Shoe consoles



Arch

The product helps evenly spread the load of your weight and well support the curved inner area of the feet, and also helps correct bent legs. It could mitigate fatigue incurred while standing or walking for a long time.



Ventilation hole

Ventilation hole helps keep you feel refreshed.



Aromatic essence

The strawberry aromatic essence is well absorbed in the product, which helps get rid of the acidic sweat while wearing shoes, and also eliminate bad smells out of shoes.




military boots



running shoes



shoes

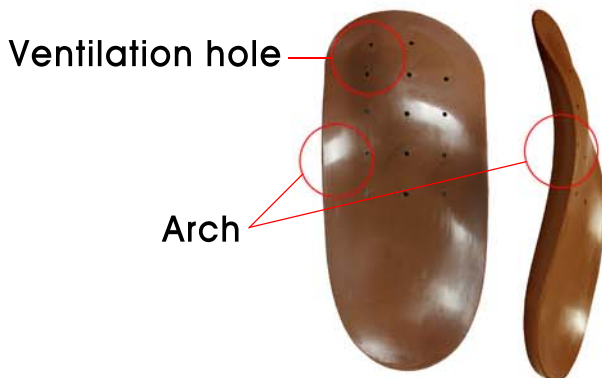

Material

P.P copolymer ,
strawberry aromatic essence



Health Bar


New Concept of arms & legs exercise equipment

By the user simply holding the acupuncture handle and turning the bar, it helps to exercise wrist joints, elbow joints, and shoulder joints smoothly and improves one's blood circulation. Also, rolling the acupuncture handle under the soles of the user's feet on the ground helps sole acupuncture and exercise of the knee



The body heat distribution using the Infrared body heat photographing


The body heat distribution of the palm of a hand.

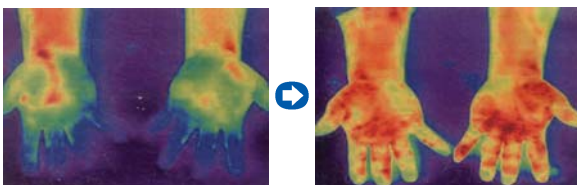

Before use

After use


The body heat distribution of the sole of a foot.

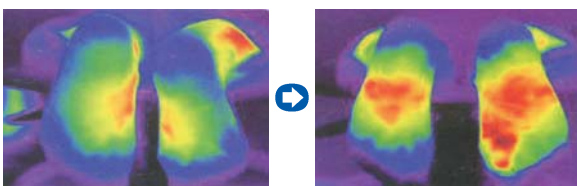


Before use

After use

Product main component

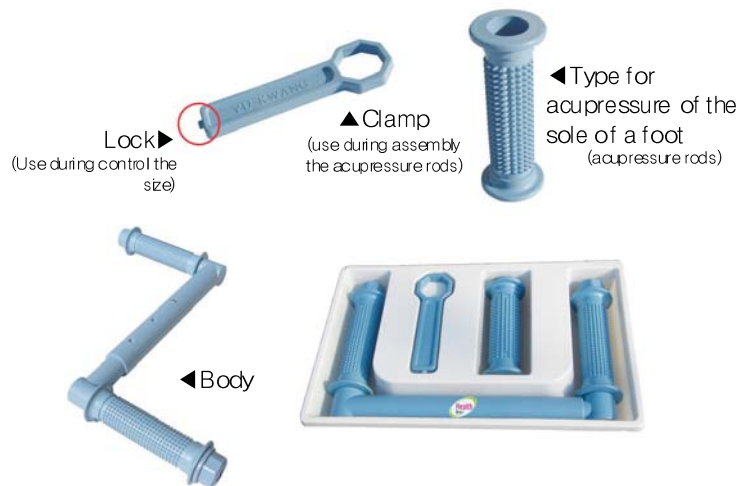

How to use example

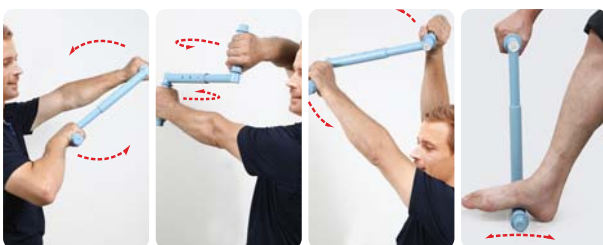

Elbow joint

Shoulder joint

Wrist joint

The sole of a foot acupuncture

Product specifications


Size	370(W) × 175(L) × 180(H)mm
Weight	840g
Material	ABS copolymer
Color	blue

Tourmaline Scalp, Skin Massager

Tourmaline is known to generate
A microcurrent of 0.06mA which is similar to that of human neurotransmissions that penetrate acupuncture points and stimulate body cells

4–11 micron far infrared ray believed to be responsible for growth

Negative Ion

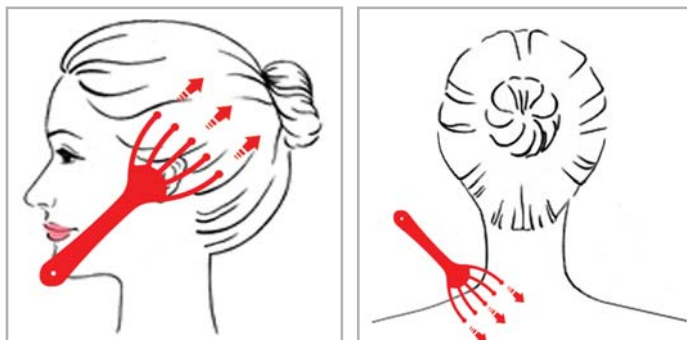

Scrub naturally on the scalp, face, back, shoulder, neck, knee, etc.

**Precautions: Do not apply too much force.*

Product specifications

Size	100(W) × 250(L) × 50(H)mm
Weight	50g
Material	P,P copolymer, Tourmaline



Scalp Masager

Skin Masager

Knee Immobiliser with Detachable Patella Strap



- It is necessary to measure the correct length of the knee immobiliser before the actual application. A correctly sized immobilizer fits on the leg 2" to 3" (5.08cm–7.62cm) above the lateral malleolus and 6" to 8" (15.24cm–20.32cm) below the femoral head.

How to assemble



- 1 Open and position "Knee Immobiliser with Detachable Patella Strap" on the leg so that the contoured posterior strays are positioned correctly.



- 2 Secure alternating loop and lock straps.



- 3 Tighten them for a secure fit.



- 4 Attach optional patella strap to help prevent swelling and to hold dressing in place. Secure patella strap over the patella

Product specifications

(Size)	width	height
소	600mm	43mm
중	600mm	50mm
대	600mm	55mm
재질	Polyester, Aluminium	



Welfare medical device



senior-friendly products



Korea Certification

Stick YGM04

This adjustable cane is helpful for seniors wanting to spend a day outside without the help of an assistant




Precautions


- Adjust length before using, according to the user's body
- Check every part before using, If rubber tip on the bottom of the cane is worn out, replace before use
- Be extra careful when using near slippery areas
- Storage : avoid hot temperatures
- People who exceed the safe operation load(90kg) should not use



Holding the stick

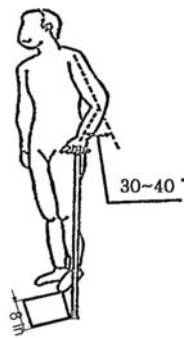

The length of the right stick

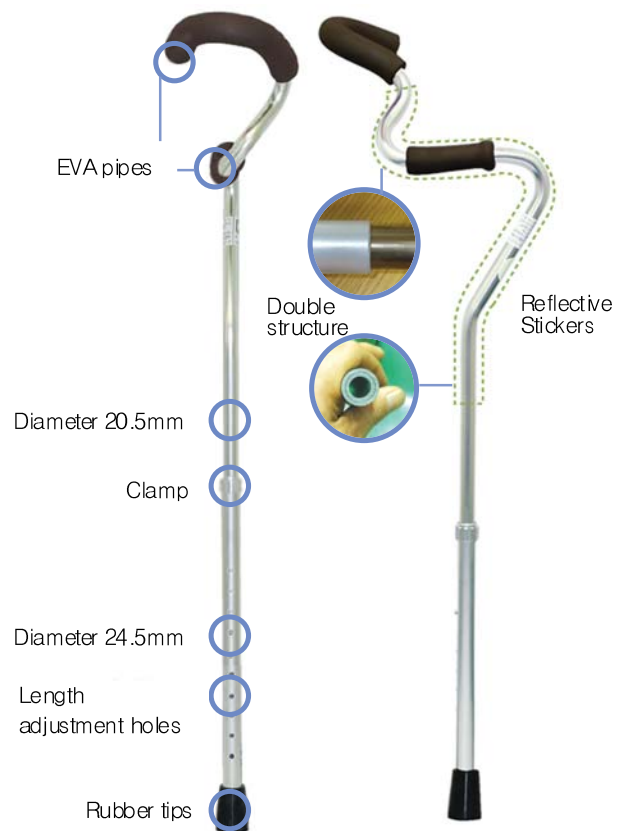

Product specifications

Size	Maximum 110cm / Minimum 88cm
Weight	800g
Material	Aluminum rods , Rubber tips, EVA pipes, Reflective stickers
Capacity	90kg





Weiba medical device



senior-friendly products



Group standard



Toilet/Shower room Grab bar

YGM03



To prevent slipping in the bathroom

This product aids seniors or people with disabilities in the shower or bathroom. It supports the users and prevents accidents, and it is simple and convenient to use. At a typical bathroom, where shower and toilet share the same space, using YGW-L and YGW-T together can be more effective.



Product specifications

Size	470(W) × 410(L) × 310(H)mm—Minimum
	500(W) × 410(L) × 310(H)mm—Maximum
Weight	2,680g
Capacity	90kg
Material	Stainless, Silicone, Aluminum

Apartment Balcony & Front Entrance/Hospital Room



• Door knob type • YGA001

Size	360mm/φ32(Bar diameter)
Weight	350g
Component	Grab Bar , piece (NICA tot 6mm 6place, piece 4mm6dace)
Material	Stainless, ABS copolymer
Capacity	100kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



• Door knob type • YGA002

Size	460mm/φ32(Bar diameter)
Weight	430g
Component	Grab Bar , piece (NICA tot 6mm 6place, piece 4mm6dace)
Material	Stainless, ABS copolymer
Capacity	100kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



• Door knob type • YGA003

Size	660mm/φ32(Bar diameter)
Weight	600g
Component	Grab Bar , piece (NICA tot 6mm 6place, piece 4mm6dace)
Material	Stainless, ABS copolymer
Capacity	100kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



• Door knob type • YGA003

Size	660mm/φ32(Bar diameter)
Weight	600g
Component	Grab Bar , piece (NICA tot 6mm 6place, piece 4mm6dace)
Material	Stainless, ABS copolymer
Capacity	100kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

**Glow-in-the-dark
Stickers on every product**



**Apartment balcony, front entrance,
glass door, and hospital sliding door,
because of their functionality...**

This product can be installed anywhere, in any arrangement—on hospital sliding door, apartment balcony, and staircase or front porch of a house, in horizontal or vertical arrangement. This Grab Bar is made especially for seniors and children, to provide convenience and support while standing or moving to another place.

Apartment Balcony & Front Entrance/Hospital Room



• Door knob type • YGH02

Size	830mm/ ϕ 32(Bar diameter)
Weight	760g
Component	Grab Bar, piece(NICA bot 6mm ϕ base, piece 4mm ϕ base)
Material	Stainless, ABS copolymer
Capacity	100kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



• Door knob type • YGH002

Size	830mm/ ϕ 32(Bar diameter)
Weight	760g
Component	Grab Bar, piece(NICA bot 6mm ϕ base, piece 4mm ϕ base)
Material	Stainless, ABS copolymer
Capacity	100kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



• Door knob type • YGA02

Size	460mm/ ϕ 32(Bar diameter)
Weight	450g
Component	Grab Bar, piece(NICA bot 6mm ϕ base, piece 4mm ϕ base)
Material	Stainless, ABS copolymer
Capacity	100kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



• Door knob type • YGC01

Size	550mm/ ϕ 32(Bar diameter)
Weight	690g
Component	Grab Bar, piece(NICA bot 8mm ϕ base, piece 4mm ϕ base)
Material	Stainless, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

How to use example



For Not only safely, but also interior aesthetics

• Yu Gwang Precision Grab Bar is designed to harmonize with any interior spaces, by upgrading existing straight Grab bars to different curves and color. Function of this Grab Bar includes: a close contact between the hand and the Grab Bar, a non-slip grip and stability. This product can be applied to a variety of uses, therefore it is worth installing.


Toilet/Shower room



• Shower & toilet type • YGC001

Size	500mm/Φ 32(Bar diameter)
Weight	640g
Component	Grab Bar, piece(NICA bot 8mm 6place, piece 4mm 8place)
Material	Stainless, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

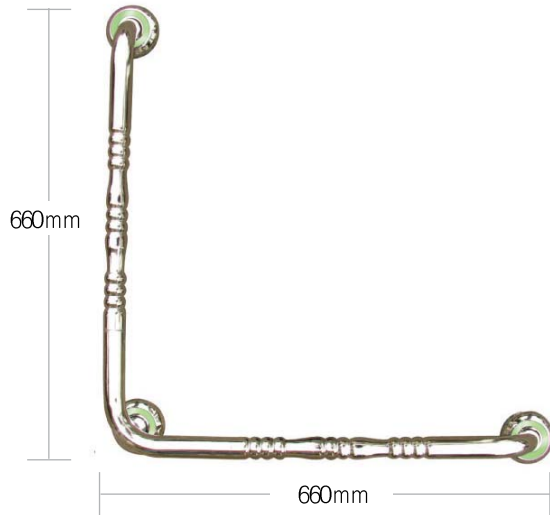

• Shower & toilet type • YGW-L

Size	660mm * 660mm/Φ 32(Bar diameter)
Weight	1500g
Component	Grab Bar, piece(NICA bot 8mm 9place, piece 4mm 12place)
Material	Stainless, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
 When attaching the plywood, please use 5mm pieces in it of over 10mm thickness immediately.

How to use example



To prevent slipping in the bathroom

- This product aids seniors or people with disabilities in the shower or bathroom. It supports the users and prevents accidents, and it is simple and convenient to use. At a typical bathroom, where shower and toilet share the same space, using YGW-L and YGW-T together can be more effective.



Welfare medical device



Senior-friendly products



Group standard

Toilet/Shower room



• Shower & toilet type • YGW-T

Size	570mm * 700mm / ϕ 32(Bar diameter)
Weight	1400g
Component	Grab Bar, piece(NICA bot 8mm 6place, piece 4mm 8place)
Material	Stainless, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep.
(drill size: 8mm)

• Shower & toilet type • WA-T100

Size	650mm * 700mm / ϕ 32(Bar diameter)
Weight	1750g
Component	Grab Bar, piece(NICA bot 8mm 9place, piece 4mm 12place)
Material	Stainless, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep.
(drill size: 8mm)

To prevent slipping in the bathroom

• This product aids seniors or people with disabilities in the shower or bathroom. It supports the users and prevents accidents, and it is simple and convenient to use. At a typical bathroom, where shower and toilet share the same space, using YGW-L and YGW-T together can be more effective.





Weiba medical device




senior-friendly products



Group standard


Silicon Grab Bar

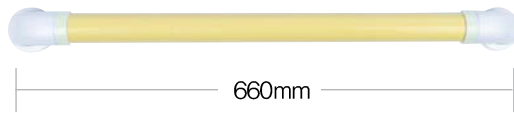

- Door knob type
- HA01

Size	660mm/Φ36(Bar diameter)
Weight	725g
Component	Grab Bar , piece(NICA bot 6mm 6place, piece 4mm 6place)
Material	Stainless, ABS copolymer, silicon
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

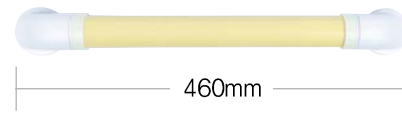

- Door knob type
- HA02

Size	460mm/Φ36(Bar diameter)
Weight	510g
Component	Grab Bar , piece(NICA bot 6mm 6place, piece 4mm 6place)
Material	Stainless, ABS copolymer, silicon
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



- Door knob type
- HA03

Size	660mm* 660mm/Φ36(Bar diameter)
Weight	1,550g
Component	Grab Bar , piece(NICA bot 8mm 6place, piece 4mm 12place)
Material	Stainless, ABS copolymer, silicon
Capacity	80kg


When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

For Not only safety, but also interior aesthetics.

•Yu Gwang Precision Grab Bar is designed to harmonize with any interior spaces, by upgrading existing straight grab bars to different curves and color. Function of this Support Bar includes: a close contact between the hand and the grab Bar, a non-slip grip and stability. This product can be applied to a variety of uses, therefore it is worth installing.




Senior housing

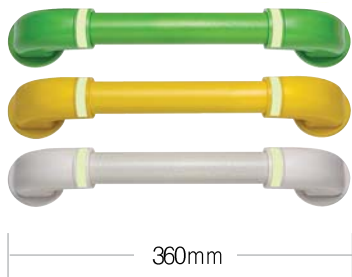

- Hospital & Senior housing facilities type
- AB001

Size	360mm/φ32(Bar diameter)
Weight	390g
Component	Grab Bar , piece (NICA bot 6mm 6/face, piece 4mm 6/face)
Material	Steel, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

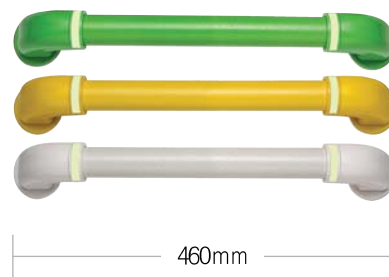

- Hospital & Senior housing facilities type
- AB002

Size	460mm/φ32(Bar diameter)
Weight	510g
Component	Grab Bar , piece (NICA bot 6mm 6/face, piece 4mm 6/face)
Material	Steel, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.


High Quality Grab Bar

- This product improves a safety issue by modifying the typical straight Grab bar to a non-slip “W” shape. By solving the aesthetic issue, the variety of color choices and different orientations of the product blends well with any surrounding environment. Also, every product is powder coated with semi-hammonton to provide the luxurious look and the feel, compared to existing products.



◀Semi-Hammonton Power Coating.

Non-Slip feature

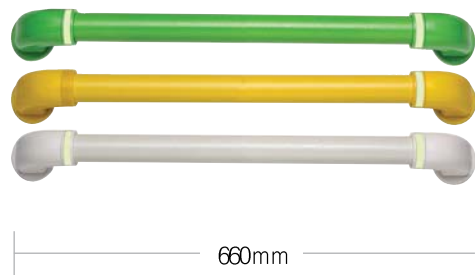

- Hospital & Senior housing facilities type
- AB003

Size	660mm/φ32(Bar diameter)
Weight	730g
Component	Grab Bar , piece (NICA bot 6mm 6/face, piece 4mm 6/face)
Material	Steel, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

Hospital



460mm

- Hospital & Senior housing facilities type
- AB2

Size	460mm/φ32(Bar diameter)
Weight	510g
Component	Grab Bar, piece(NICA bot 6mm 6place, piece 4mm 6place)
Material	Steel, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.



660mm

- Hospital & Senior housing facilities type
- AB3

Size	660mm/φ32(Bar diameter)
Weight	730g
Component	Grab Bar, piece(NICA bot 6mm 6place, piece 4mm 6place)
Material	Steel, ABS copolymer
Capacity	80kg


When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

Hospital, because of its functionality...

- At hospitals, health and safety of patients is a priority, and preventing them from falling is very important. This product is used to support a patient when going to bed, getting up, and going to a bathroom. Installing this Grab bar in each patient's room, bathroom, and hallway will maximize safety, as well as provide an aesthetically pleasing environment.




Senior housing

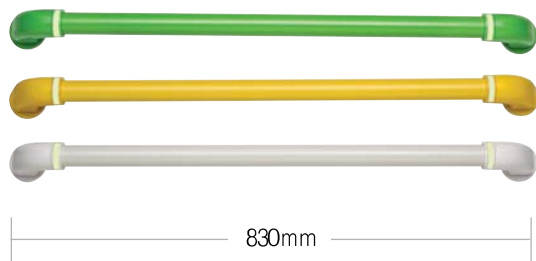

- Hospital & Senior housing facilities type
- HB002

Size	830mm/φ32(Bar diameter)
Weight	930g
Component	Grab Bar, piece(NICA bot 6mm 6place, piece 4mm 6place)
Material	Steel, ABS copolymer
Capacity	80kg


When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

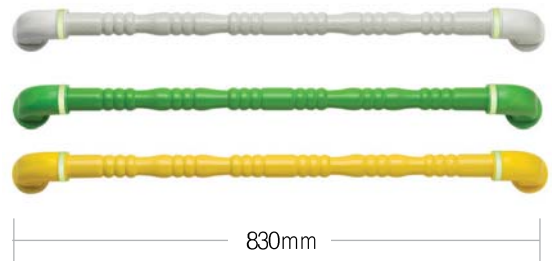

- Hospital & Senior housing facilities type
- HB2

Size	830mm/φ32(Bar diameter)
Weight	900g
Component	Grab Bar, piece(NICA bot 6mm 6place, piece 4mm 6place)
Material	Steel, ABS copolymer
Capacity	80kg


When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 6mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

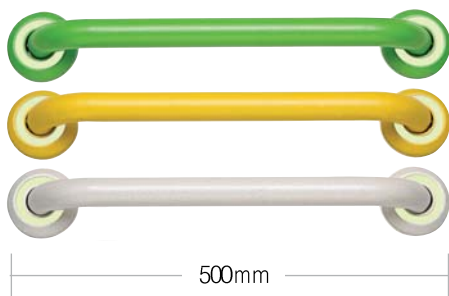

- Factory & Toilet type
- CA001

Size	500mm/φ32(Bar diameter)
Weight	710g
Component	Grab Bar, piece(NICA bot 8mm 6place, piece 4mm 6place)
Material	Steel, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

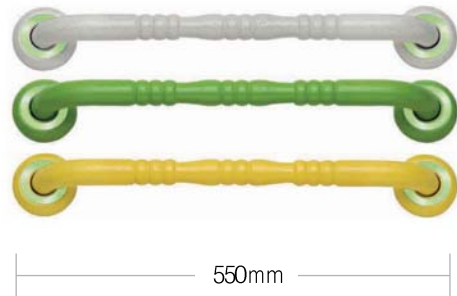

- Factory & Toilet type
- CA1

Size	500mm/φ32(Bar diameter)
Weight	740g
Component	Grab Bar, piece(NICA bot 8mm 6place, piece 4mm 6place)
Material	Steel, ABS copolymer
Capacity	80kg

When attaching to a concrete wall, drill holes that are 45mm deep. (drill size: 8mm)
 When attaching the plywood, please use 4mm pieces in it of over 10mm thickness immediately.

Innovative seedling device Porter Master

P O R T E R M A S T E R



Making hole for seedling and soil covering at once!

This Porter Master has been developed to solve the conventional problems (Difficulty in making constant holes for seedling and low productivity) and it makes seedling holes on ridge very quickly and regularly as well as recovering up with soil easily and efficiently. Further this product can be utilized in not only seedling but sowing, pot planting and weekend farm.



Sowing



Planting Seedling



Weekend Farm



Outing Class

Instruction



1 Press with thumb to spread tip.



2 Press onto ridge to make seedling hole.



3 Insert seedling into hole



4 Drop soil to cover.

Notes

- * Avoid to use on hard soil, or it may cause damage resulting in injury.
- * Child shall uses this product under adult person's supervision to prevent any accident.

Mole Spin Porter master

(Semi-automatic stick excavator)

Tough and arduous digging!
Make it easy and fast by using an electric drill !!
When sowing potatoes, peppers, corn, melons, watermelons,
and more, use the mole spin-porter Master for your convenience!

※ Compatible all type of electric drill



▲ Example of general soil

▲ Example vinyl

■ How to use

1. Mount the mole top shaft into the electric drill and fix it tightly.
2. Place a Porter master on the ground and operate the drill to drill to the desired depth.
3. Clean it with water after use.
(Stainless material 430)

■ Precautions for use

- ※ Make sure to wear gloves when working.
- ※ Fix the porter master and the electric drill tightly.
- ※ It may be sharp depending on the speed of the drill, so use it carefully.
- ※ Keep out of reach of infants and children.
- ※ The manufacturer cannot take responsibility for injuries caused by carelessness during work.



▲ Electric drill is not included



(Semi-automatic stick excavator)

Mole Spin Porter master II

Tough and arduous digging!

Make it easy and fast by using an electric drill !!

When sowing potatoes, peppers, corn, melons, watermelons, and more, use the mole spin-porter Master for your convenience!

※ Compatible all type of electric drill

Electric drill is not included ▶



▲ Example of general soil



▲ Example vinyl

■ How to use

1. Adjust the handle height to fit your height, and adjust spacer the height and width of the gap to the desired width.
2. Mount the mole top shaft into the electric drill and fix it tightly.
3. Place a Porter master II on the ground and operate the drill to drill to the desired depth.
4. Clean it with water after use. (Stainless material 430)

■ Precautions to use

- ※ Make sure to wear gloves when working
- ※ Fix the porter master and the electric drill tightly.
- ※ It may be sharp depending on the speed of the drill, so use it carefully.
- ※ Keep out of reach of infants and children.
- ※ The manufacturer cannot take responsibility for injuries caused by carelessness during work.





Venture Company
Yu Gwang Precision

Venture Company
Yu Gwang Precision

85, Gongdan Ro 3, Buk-gu, Daegu, KOREA 41497
Tel. +82-53-357-5512~3- Fax. +82-53-353-5543
E-mail. injo4162@naver.com, info@ukm.co.kr Homepage : www.ukm.co.kr